



PROMOTIONS | 2026
DU 1^{ER} SEPTEMBRE AU 31 DÉCEMBRE

Offres de fin d'année
AUTOMNE **HIVER**

L'EXCELLENCE AU SERVICE DES CABINETS DENTAIRES



Matériel premium, design moderne, décliné en plusieurs couleurs pour s'adapter à votre environnement. Avec Prodentalpes, offrez à votre cabinet des équipements alliant performance, confort et esthétique.

NOS ATOUTS

- ✓ **Qualité et fiabilité** – Des produits sélectionnés pour leur durabilité et leur efficacité
- ✓ **Confort et ergonomie** – Pensés pour améliorer l'expérience des praticiens et des patients.
- ✓ **Service expert et réactif** – Un accompagnement personnalisé pour répondre à vos besoins.

CONTACTEZ-NOUS

✉ contact@prodentalpes.fr

☎ 06 01 21 16 35



prodentalpes.fr

ZI DES DRAGIEZ, 30 RUE DE GRÉBELIN,
74800 LA ROCHE-SUR-FORON


OFFRE INSERTS 2026

PROMOTION 10+2 OFFERTS

✓ AVANTAGES :

Vous pouvez retrouver tous les inserts Satelec, EMS, Mectron, Kavo, NSK, Sirona et Amdent, sauf pour les piezosurgery.
(Mectron 90€ au lieu de 160€)

DE 39€ À 69€

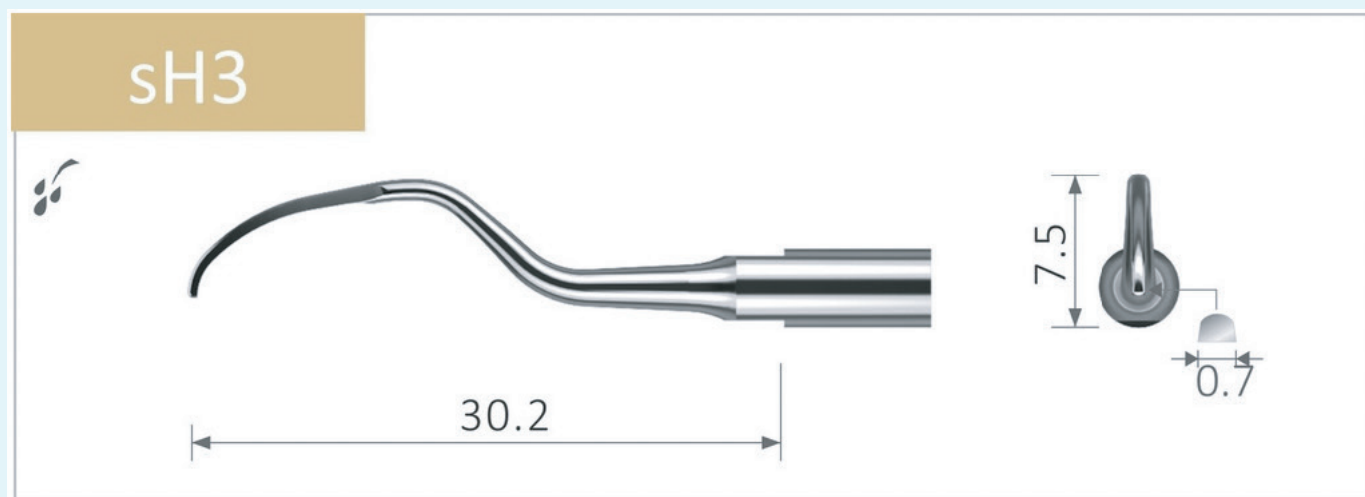

✓ PRIX AVANTAGE PRODENTALPES

Promotion : H3 à 39 €
Promotion : IS à 39 €

FABRICANT AUTRICHIEN 



SCALING

Scaling is precise work. Proper instrument selection is essential to achieving complete periodontal cleaning. Our wide variety of tips give you the freedom to choose the exact instrument for the situation.



PERIODONTICS

Perio tips are thin and designed for root planing and maintenance to provide the best access to furcation and curved roots.



ENDODONTICS

Endodontic tips can be used in many areas in endodontics. They are excellent for the removal of posts, removing dentin in pulp chambers, finding and widening orifices, preparing canals, removing broken instruments and cleaning prepared canals.



CAVITY PREPARATION

These tips are diamond coated and can be used to prepare a cavity in the tooth before carrying out further dental work.



SURGERY TIPS

Surgery tips are designed for a wide range of bone surgery applications, including sinus lifting and implantation procedures.



Principle of Operation.

For piezo style tips, including all Xpedent tips except those to use with NSK® and Kavo Soniflex® compatible air scalers, power is generated by piezoelectric crystals. For airscalers power is generated by compressed air. This creates a sinusoidal motion with the greatest power (amplitude) at the end of the tip followed by a weaker area, then a stronger area in accordance with a sine wave. The effect of the tip is therefore very dependent upon the part of the tip in contact with the tooth. Contact with the side of the tip in the second strong power band, 2-3 mm from the end, creates the most effective action.

The ideal power setting will depend on the tip used and the condition of the area to be treated. Endo and Perio work requires longer and thinner tips than scaling, so lower powers are appropriate. Always start with a low power and increase gradually to maximise effectiveness and minimise patient discomfort.

Orientation and Application.

To create a combination of patient comfort and efficient scaling it is important to use the side of the tip at the second power band and at a very shallow angle to the tooth. Too great an angle will bring the point of the tip in contact with the tooth and lead to patient discomfort.


Unlike hand scaling, ultrasonic scaling requires only a very gentle touch to the tooth surface; the ultrasonics provide all the power so the dental practitioner only needs to guide the tip very gently around the area to be treated. The removal of hard deposits can be expedited by applying the end of the tip where the vibrations are strongest; this will cause calculus to break away from the tooth. Once this is done, resume gentle motions. Too much force, or over-use of the tip end, will cause discomfort.

CONTENT

Generic tips		6
Xpudent tips (suitable for use with EMS® scaler machines). The Xpudent brand belongs to a third party which has no connection with the original equipment manufacturer or its representatives.		12
Xpudent tips (suitable for use with SATELEC® scaler machines). The Xpudent brand belongs to a third party which has no connection with the original equipment manufacturer or its representatives.		16
Xpudent tips (suitable for use with NSK® scaler machines). The Xpudent brand belongs to a third party which has no connection with the original equipment manufacturer or its representatives.		21
Xpudent tips (suitable for use with SIRONA® scaler machines). The Xpudent brand belongs to a third party which has no connection with the original equipment manufacturer or its representatives.		24
Xpudent tips (suitable for use with KAVO® scaler machines). The Xpudent brand belongs to a third party which has no connection with the original equipment manufacturer or its representatives.		26
Xpudent tips (suitable for use with AMDENT® scaler machines). The Xpudent brand belongs to a third party which has no connection with the original equipment manufacturer or its representatives.		32
Xpudent tips (suitable for use with MECTRON® scaler machines). The Xpudent brand belongs to a third party which has no connection with the original equipment manufacturer or its representatives.		34

SIRONA® is a registered trademark of Sirona Dental Systems GmbH
EMS® is a registered trademark of Electro Medical Systems SA
SATELEC® is a registered trademark of ACTEON Group
Mectron® is a registered trademark of Mectron s.p.a.
Amdent® is a registered trademark of LM Instruments Oy
Kavo® is a registered trademark of KAVO Dental Corporation
NSK® is a registered trademark of Nakanishi INC. Corporation

CONTENT

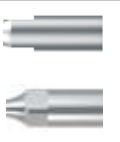Xpedent tips (suitable for use with SATELEC® surgery machines). The Xpedent brand belongs to a third party which has no connection with the original equipment manufacturer or its representatives.		41
Xpedent tips (suitable for use with EMS® surgery machines). The Xpedent brand belongs to a third party which has no connection with the original equipment manufacturer or its representatives.		41
Xpedent torque wrenches & wrenches		45
U-Files		46
Thread identification		46
Xpedent scaler handpieces		47
Autoclavable tip holders		47
Xpedent built-in scaler machines		48
Tip wallets		49

SIRONA® is a registered trademark of Sirona Dental Systems GmbH
EMS® is a registered trademark of Electro Medical Systems SA
SATELEC® is a registered trademark of ACTEON Group
Mectron® is a registered trademark of Mectron s.p.a.
Amdent® is a registered trademark of LM Instruments Oy
Kavo® is a registered trademark of KAVO Dental Corporation
NSK® is a registered trademark of Nakanishi INC. Corporation




diamond coated



irrigation

Xpedent tips suitable for use with EMS® scaling devices
 Xpedent tips suitable for use with SIRONA® scaling devices
 Xpedent tips suitable for use with SATELEC® scaling devices
 Xpedent tips suitable for use with NSK® scaling devices
 Xpedent tips suitable for use with MECTRON® scaling devices

SCALING

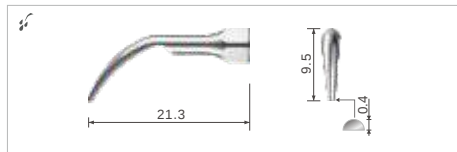

Used to remove light and medium supragingival calculus and bacterial plaque.


G1

GS1

GN1

GD1

GM1

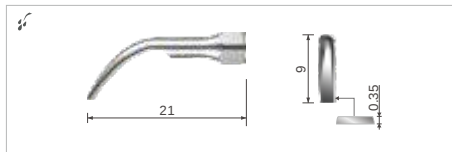


Flat edge is used to remove heavy supragingival calculus.

G2

GS2

GN2

GD2

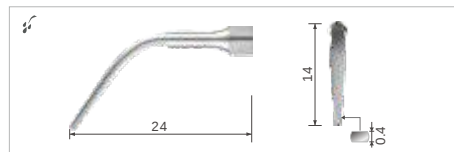

Used to remove calculus and bacterial plaque from supragingival interdental areas.


G3

GS3

GN3

GD3

GM3

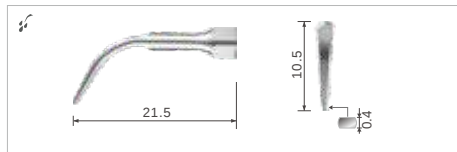

Used to remove all supragingival calculus and bacterial plaque.


G4

GS4

GN4

GD4

GM4

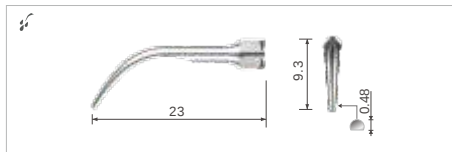

Used to remove all supragingival calculus and bacterial plaque.


G5

GS5

GN5

GD5

GM5

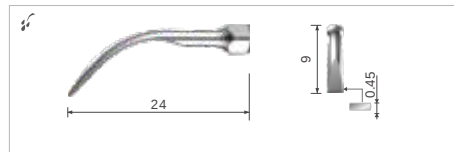

Used to remove heavy supragingival calculus.


G6

GN6

GD6

GM6

PERIO

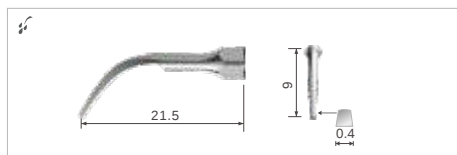

Used to remove subgingival calculus.


P1

PS1

PN1

PD1

PM1

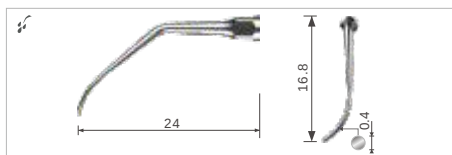


Curved to the left, for treating furcations, interdental calculus and very narrow spaces.

P2L

PN2L

PD2L

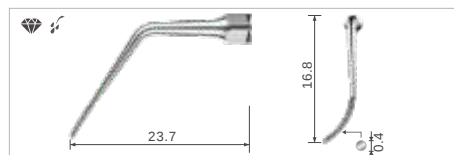


Curved to the left and diamond coated, for treating furcations, interdental calculus and very narrow spaces.

P2LD

PD2LD

SIRONA® is a registered trademark of Sirona Dental Systems GmbH
 EMS® is a registered trademark of Electro Medical Systems SA
 SATELEC® is a registered trademark of ACTEON Group
 Mectron® is a registered trademark of Mectron s.p.a.
 Amdent® is a registered trademark of LM Instruments Oy
 Kavo® is a registered trademark of KAVO Dental Corporation
 NSK® is a registered trademark of Nakanishi INC. Corporation

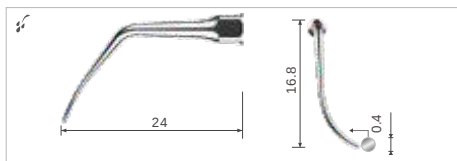


Curved to the right, for treating furcations, interdental calculus and very narrow spaces.

P2R

PN2R


PD2R

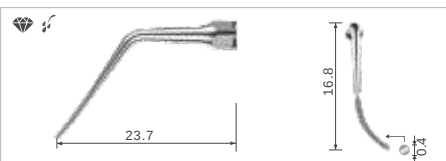

Curved to the right and diamond coated, for treating furcations, interdental calculus and very narrow spaces.

P2RD

PD2RD

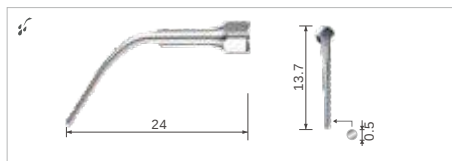

Used to remove subgingival calculus.


P3

PS3

PN3

PD3

PM3

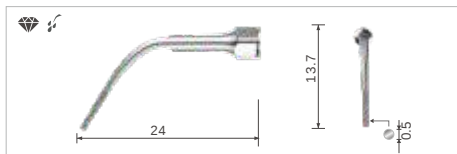


Diamond coated, used to smooth the surface of root during periodontal flap surgery.

P3D

PS3D

PN3D

PD3D

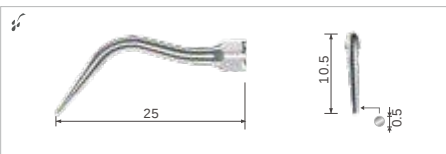


Used to remove subgingival calculus.

P4

PS4

PN4

PD4

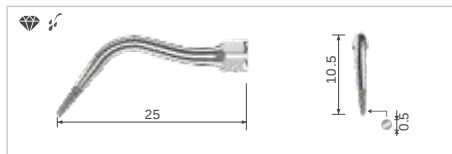


Diamond coated, a powerful tool for the removal of calcification, this tip can be used for widening calcified roots in the coronal third.

P4D

PN4D


PD4D

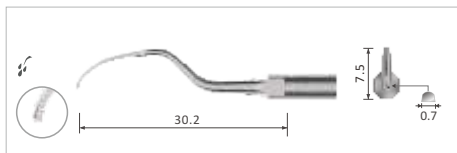

Anterior Periodontics. Enables treatment of the incisor-canine block. The guide edge is oriented parallel to the pocket. Enables work in the periodontal pocket without injury to the ligament.

eH3 (P5)


sH3

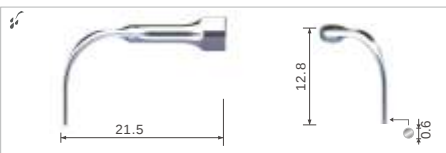

Curved to the left, used to remove subgingival calculus.

P6L

PD6L

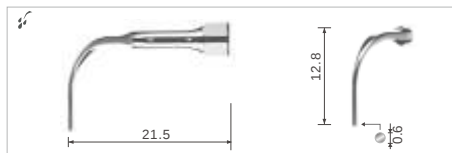

Curved to the right, used to remove subgingival calculus.

P6R

PD6R

ENDO



120° angle holder for use with Xpident U files with a shank diameter of 0.57mm, this tip is usually used for anterior root canal cleaning.

E1

ES1

EN1

ED1




95° angle holder for use with Xpident U files with a shank diameter of 0.57mm, this tip is usually used for anterior root canal cleaning.

E2

ES2

EN2

ED2

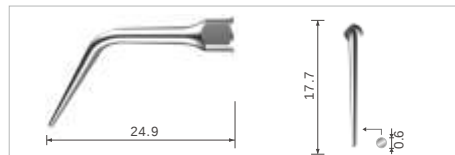


With cooling spray disabled, this tip is used for lateral thermal condensation of gutta percha.

E3

ES3

EN3

ED3

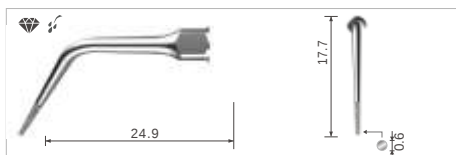


Diamond coated, used to remove calcification and bad filling material from pulp cavity and overhang pulp chamber.

E3D

ES3D

EN3D

ED3D

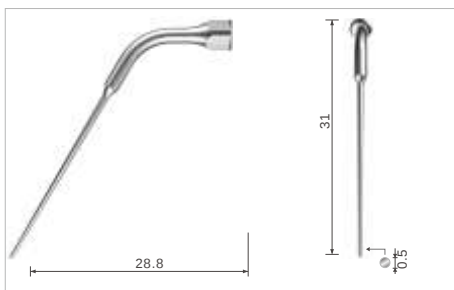


Quickly removes obstructions and broken instruments in the mid and upper thirds of large and straight canals.

E4

ES4

EN4

ED4

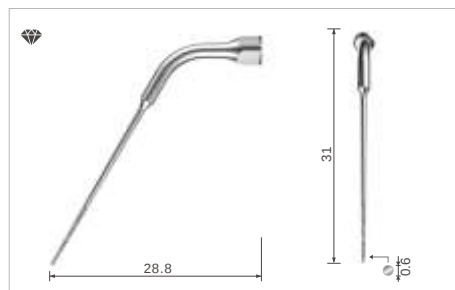


Diamond coated, for retreatments or extremely hard materials.

E4D

ES4D

EN4D

ED4D

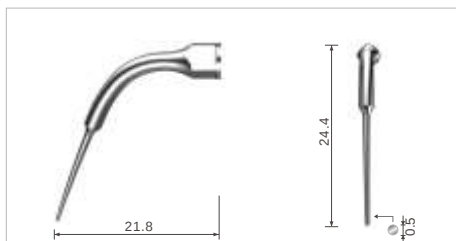


Useful for quickly eliminating calculus and temporary cements.

E5

ES5

EN5

ED5

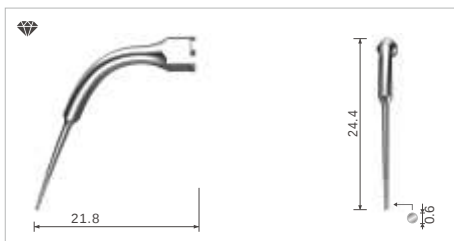

Diamond Coated, for enhanced cutting and abrasion. Especially useful for removing very hard materials on the canal wall.

E5D

ES5D

EN5D

ED5D



Quickly removes obstructions and broken instruments in the mid and upper thirds of large and straight canals.

E14

EN14

ED14




Diamond coated for retreatments or extremely hard materials.

E14D

EN14D

ED14D

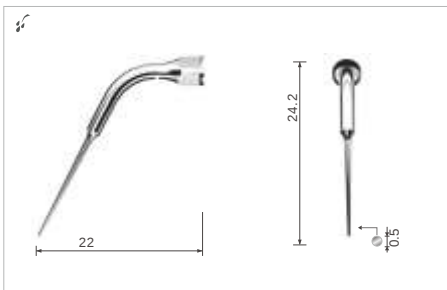


Useful for quickly eliminating calculus and temporary cements.

E15

EN15

ED15

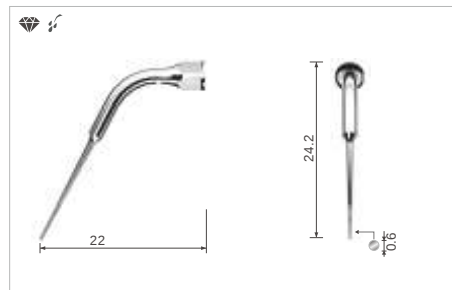


Diamond coated, for enhanced cutting and abrasion. Especially useful for removing very hard materials on the canal wall.

E15D

EN15D

ED15D

RETRO

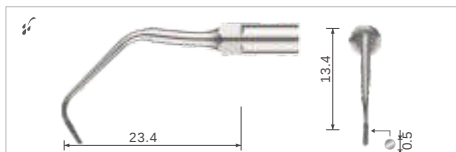


Used for gentle treatment of the root canal of anterior teeth.

E10

EN10

ED10

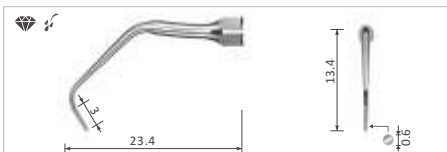


Diamond coated, used for efficient apical polishing of anterior teeth.

E10D

EN10D

ED10D

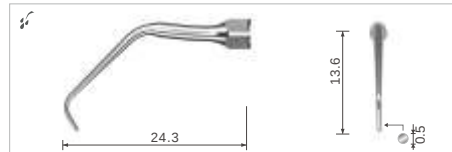


Used for gentle treatment in the root canal of anterior teeth.

E11

EN11

ED11

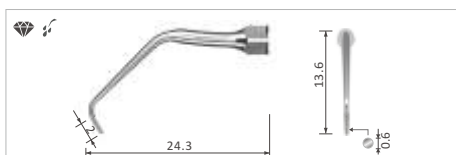


Diamond coated, with enhanced cutting and abrasive power; especially useful for removing very hard materials on the canal wall and lateral cutting.

E11D

EN11D


ED11D

CROWN REMOVAL

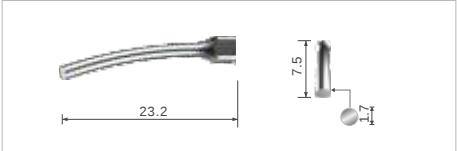

Used to remove bridges, crowns and posts by fracturing the cement.


- G7
- GS7
- GN7
- GD7

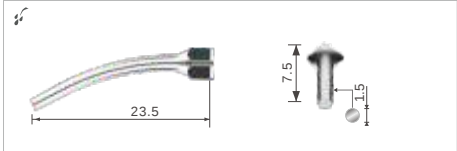

Used to remove bridges, crowns and posts by fracturing the cement.


- G8

CAVITY PREPARATION

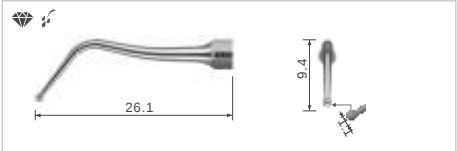


Ball diamond tip, for preparation of occlusal surfaces and cervical ridges.

- SB1
- SBS1
- SBN1
- SBD1

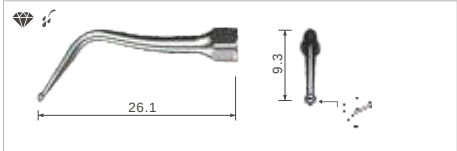


Semi-spherical diamond tip, for finishing cervical interproximal areas without causing lesions to adjacent teeth.

- SB2
- SBS2
- SBN2
- SBD2

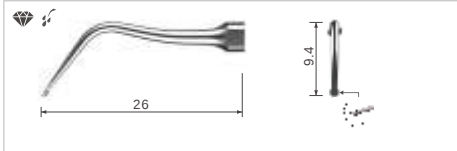

Semi-spherical diamond tip, for finishing cervical interproximal areas without causing lesions to adjacent teeth.

- SB3
- SBS3
- SBN3
- SBD3

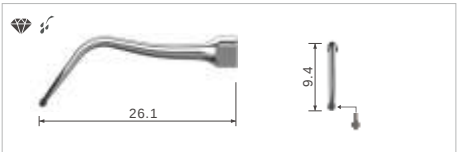

Semi-spherical diamond tip set 45° to the left. Allows the dentist to remove carious areas of a tooth without damaging the adjacent teeth.

- SBL
- SBSL
- SBNL
- SBDL



Semi-spherical diamond tip set 45° to the right. Allows the dentist to remove carious areas of a tooth without damaging the adjacent teeth.

- SBR
- SBSR
- SBNR
- SBDR

FG-BUR HOLDER



FG-bur holder used to expand the root canal and grind teeth.

E8

EN8

ED8




FG-bur holder used to expand the root canal and grind teeth.

E9

EN9

ED9

IMPLANT CLEANING

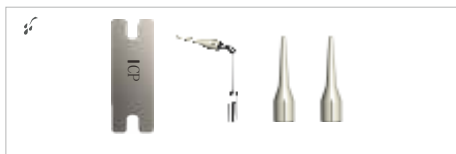

Implant Cleaning. The ICP tip is used for sub- and supragingival cleaning of implants and restorations. The replaceable peek makes the action safe and gentle. The set contains 1 ICP, 1 ICP wrench and 2 spare peeks. Replacement

eICP SET

srICP SET

nICP SET

slCP SET



Replacement. ICPpeeks are autoclavable and reusable. Replacing the ICPpeek is simple using the wrench included in ICP SET.

ICPpeek

ICPpeek

ICPpeek

ICPpeek



Xpedent tips (suitable for use with EMS® scaler machines).
The Xpedent brand belongs to a third party which has no connection with the original equipment manufacturer or its representatives.

SCALING



Universal. Accurate removal of solid supragingival calculus in all quadrants.


PERIO




Perio. The eP instrument is well suited for the precise removal of stubborn subgingival and supragingival calculus and concretions in all quadrants.



Anterior Periodontics. Enables treatment of the incisor-canine block. The guide edge is oriented parallel to the pocket. Enables work in the periodontal pocket without injury to the ligament.

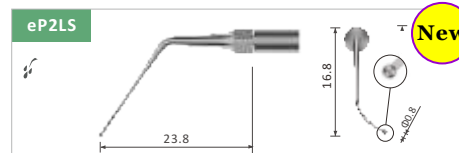


Perio. Left curved instrument with a small ball > 0.8 mm diameter on the instrument tips for periodontal debridement of furcations and concavities

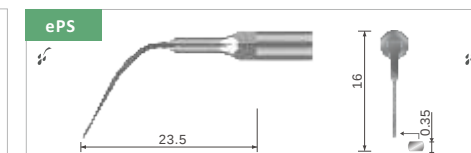


PerioSlim. The fast and minimally-invasive ePS instrument is ideal for calculus removal in deep periodontal pockets in all quadrants.




Premolar and Molar. Left-curve. First instrument for treating all the surfaces and the furcations.
Maxillary. Buccal and distal surfaces of sector 2, pivots at 13, then the buccal and mesial surfaces of sector 1.
Mandibular. Buccal and distal surfaces of sector 4, pivots at 43, then lingual and mesial surfaces of sector 3.

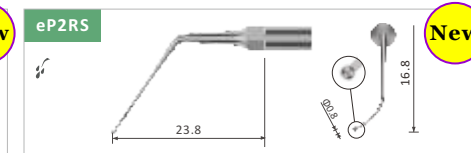


Perio. Right curved instrument with a small ball > 0.8 mm diameter on the instrument tips for periodontal debridement of furcations and concavities

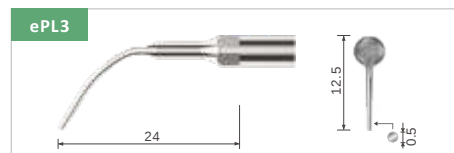

Perio Cleaning. For cleaning and disinfecting periodontal pockets, it is especially suitable for follow-up and recall visits for periodontal maintenance patients.

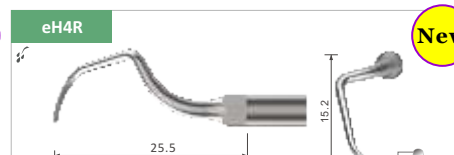


Premolar and Molar. Right-curve. Second instrument, use after sH4L.
Maxillary. Palatine and mesial surfaces of sector 2, pivots at 13, then buccal and distal surfaces of sector 1.
Mandibular. Lingual and mesial surfaces of sector 4, pivots at 43, then buccal and distal surfaces of sector 3.

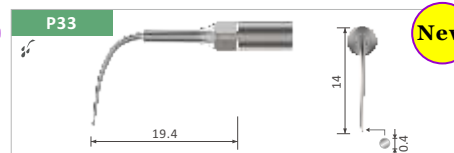


Perio. Used to remove dental calculus deep in the periodontal pocket and clean the periodontal pocket.

ENDO

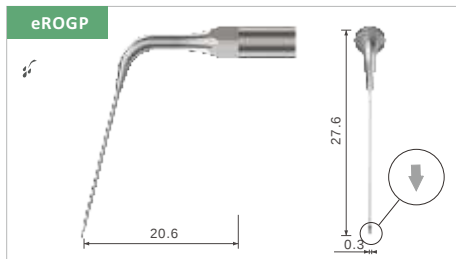

Access Long Canals. The ROGP tip is used to access areas not reached by other tips. Particularly effective for retreatments, especially in long oval canals, where most other tips would not fit. Excellent for Gutta-Percha removal.

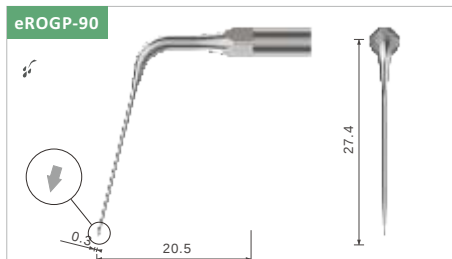

Access Long Canals - 90°. The ROGP tip is used to access areas not reached by other tips. Particularly effective for retreatments, especially in long oval canals, where most other tips would not fit. Excellent for Gutta-Percha removal.

CONDENSATION



Condensation. For cementing ceramic inlays, onlays and veneers. The eC20 SET contains 1 eC20 and 2 spare peeks.



Replacement. C20peeks are autoclavable and reusable.

IMPLANT CLEANING




Implant Cleaning. The ICP tip is used for sub- and supragingival cleaning of implants and restorations. The replaceable peek makes the action safe and gentle. The set contains 1 ICP, 1 ICP wrench and 2 spare peeks.



Replacement. ICPpeeks are autoclavable and reusable. Replacing the ICPpeek is simple using the wrench included in ICP SET.


SETS (TORQUE WRENCH INCL. SCALER TIP)

TWT eA

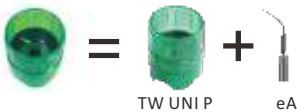


TW UNI P eA

TWT eP

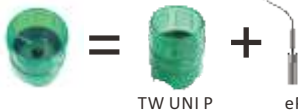

TW UNI P eP

TWT ePS

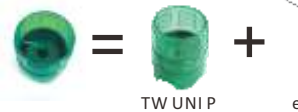

TW UNI P ePS

TWT ePL3



TW UNI P ePL3

TW UNI P



Universal
Plastic Torque Wrench




16




Xpedent tips (suitable for use with SATELEC® scaler machines).
The Xpedent brand belongs to a third party which has no connection with the original equipment manufacturer or its representatives.

SCALING

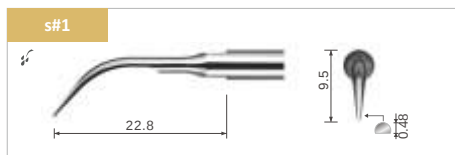


Universal. Simple supragingival scaling and tartar removal without damage to the enamel.

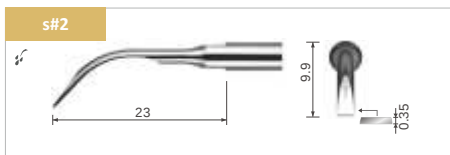


Calculus. Apply the flat part to the tooth to remove major calculus deposits.




Stain Removal. Rounded leading edge is effective for the removal of marks or stains such as tobacco or tea.

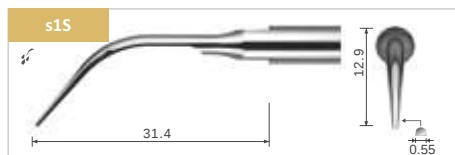


Slim. Supra and subgingival scaling. Particularly effective for interproximal spaces.

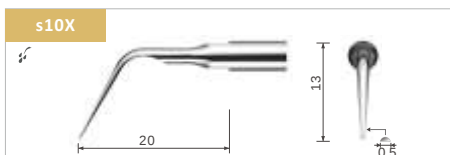

Interproximal. Its flat active part makes it suitable for the interproximal spaces and supragingival scaling. Its anatomical design enables fast and efficient use.

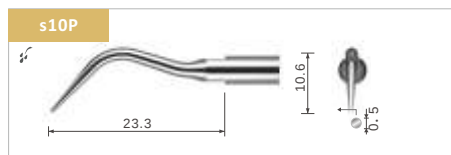

Shallow. Slim round instrument suitable for scaling shallow pockets less than 2-3mm deep.

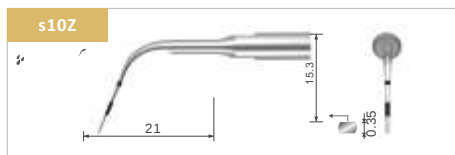


Subgingival. Medium pockets (< 4mm). Removal of biofilm and soft deposits, enables observation of pocket depth using the 3mm scale.


PERIO



Anterior Periodontics. Enables treatment of the incisor-canine block. The guide edge is oriented parallel to the pocket. Enables work in the periodontal pocket without injury to the ligament.

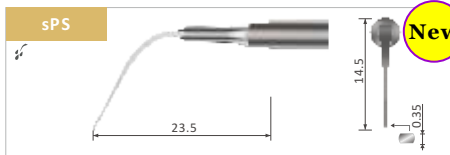

PerioSlim. The fast and minimally-invasive sPS instrument is ideal for calculus removal in deep periodontal pockets in all quadrants.

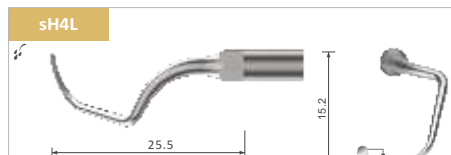


Premolar and Molar. Left-curve. First instrument for treating all the surfaces and the furcations.

Maxillary. Buccal and distal surfaces of sector 2, pivots at 13, then the buccal and mesial surfaces of sector 1.

Mandibular. Buccal and distal surfaces of sector 4, pivots at 43, then lingual and mesial surfaces of sector 3.

sH4R

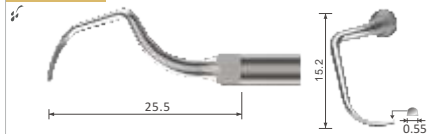


Premolar and Molar. Right-curve. Second instrument, use after sH4L.

Maxillary. Palatine and mesial surfaces of sector 2, pivots at 13, then buccal and distal surfaces of sector 1.


Mandibular. Lingual and mesial surfaces of sector 4, pivots at 43, then buccal and distal surfaces of sector 3.

sTK1L

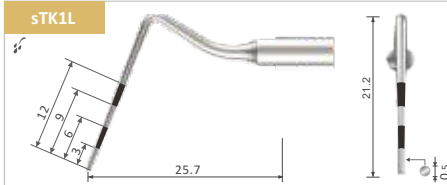

Extended Probe. Ideal for examining and maintenance of medium to deep pockets (> 4mm). A diagnosis aid used during the debridement and irrigation of pockets.

sTK1S

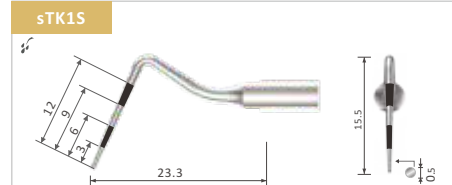


Probe. Recommended for examining shallow and medium pockets (< 4mm) and for the maintenance of simple cases. sTK1S tip is graduated every 3mm.

sTK2L



Left-Oriented Maintenance. Recommended for the maintenance of moderate to deep pockets and furcations. Equivalent to the Nabers probe.


sTK2R

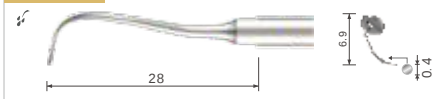

Right-Oriented Maintenance. Recommended for the maintenance of moderate to deep pockets and furcations. Equivalent to the Nabers probe.

ENDO


sET20D

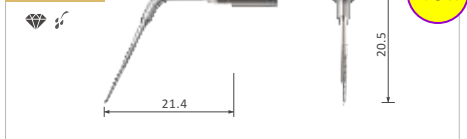

Diamond Coated Retreatment. Used for the removal of hard materials in the 1st coronal third, 30 µm diamonds,

New

sET25

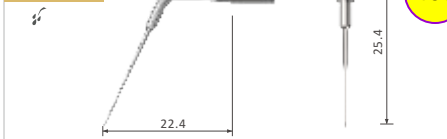

Retreatment. Preferred to use in the middle and apical thirds for retreatment and the extraction of broken instruments, Titanium-alloy, length 20mm, taper 3%.


New

sET25S




Retreatment Short. Used in coronal third and the isthmuses, Titanium-alloy, length 15mm, taper 4%.


New

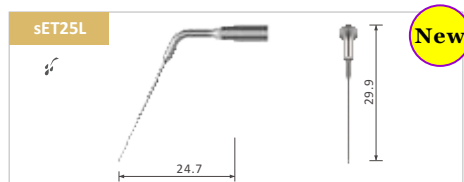


Retreatment Long. Used in the apical third and long straight canals for retreatment, Titanium-alloy, length 25mm, taper 3%.

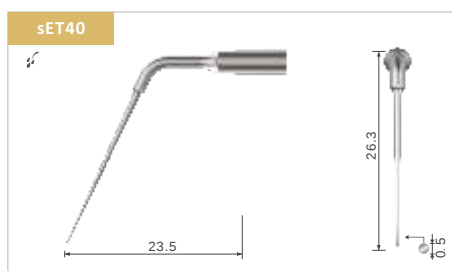


Instrument Removal. Used in the middle third of wide, straight canals, the sET40 is a rapid solution for removal of broken instruments.

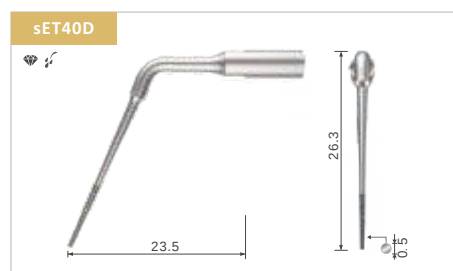


Diamond. A diamond-coated steel tip for retreatment of very hard material in the middle third. Length 40mm, 30 μ m diamonds, taper 4%.

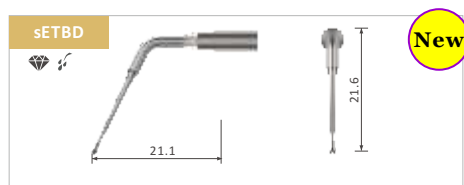

Calicification Locator. Used for searching for canals and locating calified canals, diamond coated ball, length 20mm, taper 5 %.

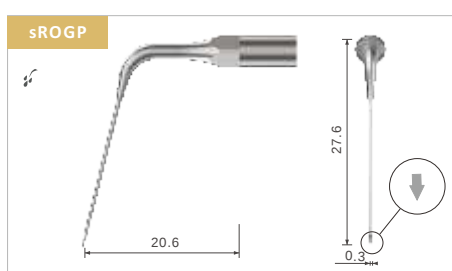


Access Long Canals. The ROGP tip is used to access areas not reached by other tips. Particularly effective for retreatments, especially in long oval canals, where most other tips would not fit. Excellent for Gutta-Percha removal.

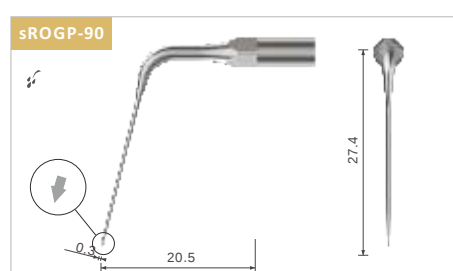


Access Long Canals - 90°. The ROGP tip is used to access areas not reached by other tips. Particularly effective for retreatments, especially in long oval canals, where most other tips would not fit. Excellent for Gutta-Percha removal.


RETRO

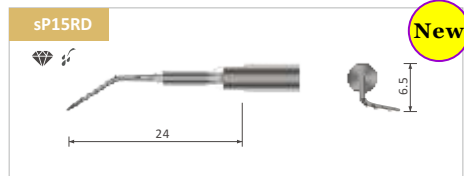


Retro Surgery Tip - Right. Used for preparation of premolars and molar canals, length 5mm, 30 μ m diamonds, taper 7%.

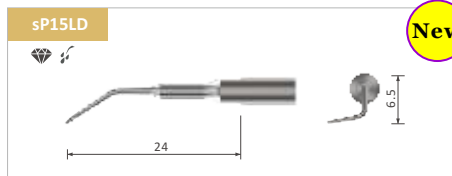


Retro Surgery Tip - Left. Used for preparation of premolars and molar canals, length 5mm, 30 μ m diamonds, taper 7%.

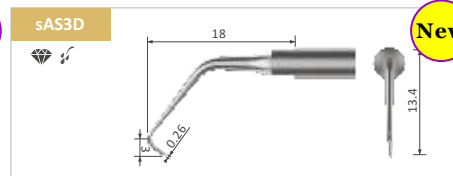

Universal Retro Surgery Tip. Used with minimum power and low pressure for apical surgery of anterior teeth, length 3mm, 30 μ m diamonds, taper 9%.

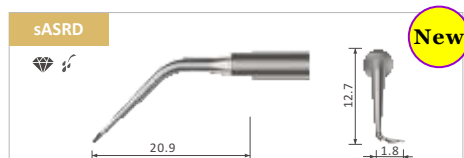

Retro Surgery Tip - Right. Used with low pressure for apical surgery of molars and premolars, length 3mm, 30 μ m diamonds, taper 10%.

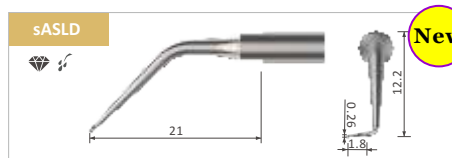

Retro Surgery Tip - Left. Used with low pressure for apical surgery of molars and premolars, length 3mm, 30 μ m diamonds, taper 10%.

CONDENSATION

sC20 SET



New

Condensation. For cementing ceramic inlays, onlays and veneers. The sC20 SET contains 1 sC20 and 2 C20peeks.

C20peek



10X

New

Replacement. C20peeks are autoclavable and reusable.

IMPLANT CLEANING

sICP SET



Implant Cleaning. The ICP tip is used for sub- and supragingival cleaning of implants and restorations. The replaceable peek makes the action safe and gentle. The set contains 1 ICP, 1 ICP wrench and 2 spare peeks.

ICPpeek



10X

Replacement. ICPpeeks are autoclavable and reusable. Replacing the ICPpeek is simple using the wrench included in ICP SET.

TW UNI P



Universal
Plastic Torque Wrench



Xpedent tips (suitable for use with NSK® scaler machines).
The Xpedent brand belongs to a third party which has no connection with the original equipment manufacturer or its representatives.

Xpedent tips suitable for use with VARIOS®

SCALING

nG1



Universal. Removal of supragingival and marginal calculus, removal of calculus from narrow interdental regions.

nG4



Scaling. Removal of supragingival and marginal calculus.

nG6




Slim. Removal of supragingival and marginal gingival calculus and removal of calculus from interproximal spaces.

nG8



Universal Slim. Removal of supra and marginal gingival calculus; particularly effective in interdental spaces.

nG16

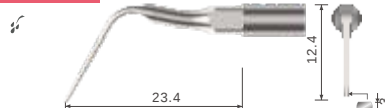

Interproximal. Removal of calculus from interproximal areas.

CONDENSATION

nC20 SET



New

Condensation. For cementing ceramic inlays, onlays and veneers. The nC20 SET contains 1 nC20 and 2 C20peaks.

C20peak



New

Replacement. C20peaks are autoclavable and reusable.

IMPLANT CLEANING

nICP SET



Implant Cleaning. The ICP tip is used for sub- and supragingival cleaning of implants and restorations. The replaceable peek makes the action safe and gentle. The set contains 1 ICP, 1 ICP wrench and 2 spare peeks.

ICPpeek



Replacement. ICPpeeks are autoclavable and reusable. Replacing the ICPpeek is simple using the wrench included in ICP SET.




TW UNI P








Universal
Plastic Torque Wrench

Xpedent tips suitable for use with TI-MAX®



SCALING

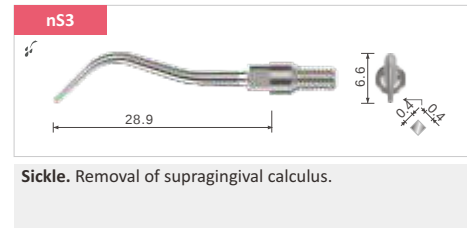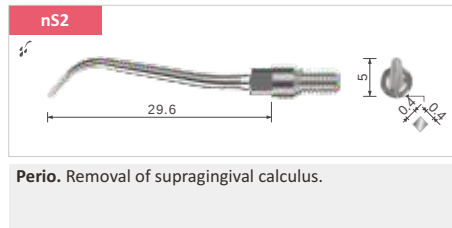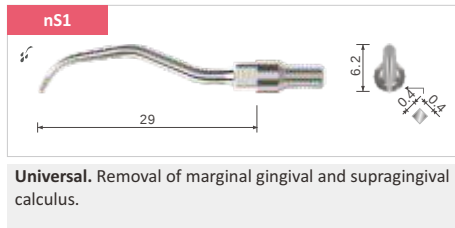

PERIO

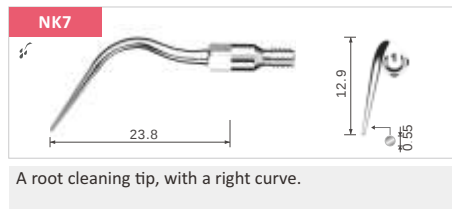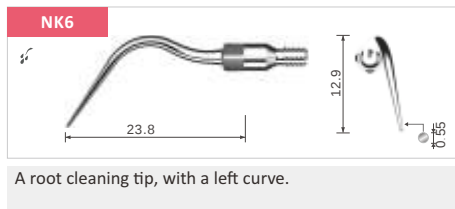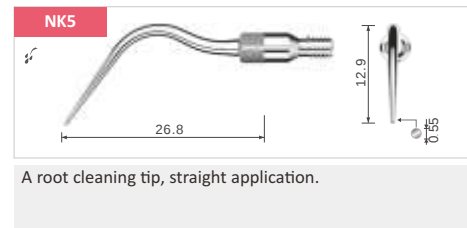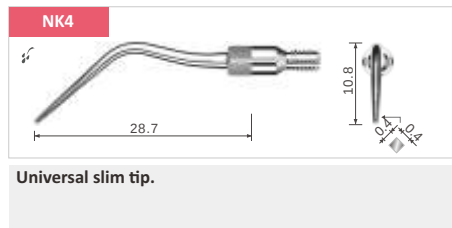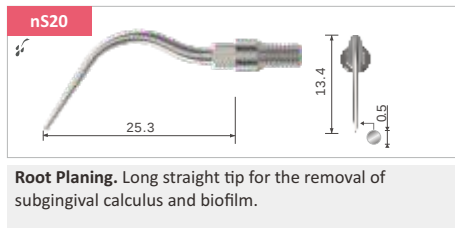


IMPLANT CLEANING

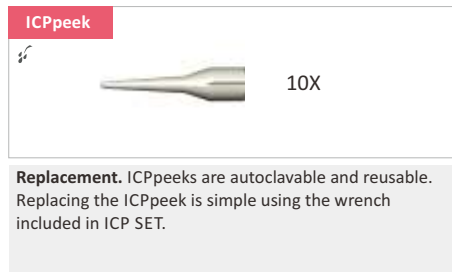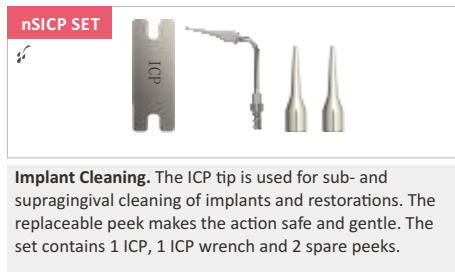




Xpedent tips (suitable for use with SIRONA® scaler machines).
The Xpedent brand belongs to a third party which has no connection with the original equipment manufacturer or its representatives.


SCALING

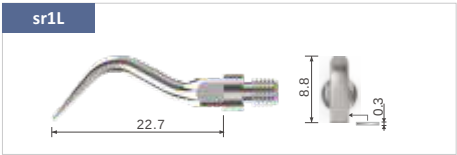

General. A flat tipped instrument for initial lingual and buccal scaling.

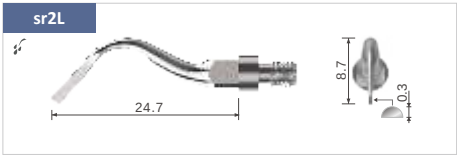

Versatile. A versatile instrument used for lingual, buccal and approximal scaling.

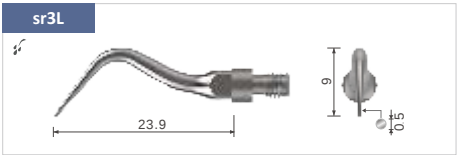


Interdental. A universal instrument which can be used for subgingival scaling at interdental sites.




Universal. A fine instrument, ideal for supragingival and subgingival tooth surfaces and interdental areas.


PERIO

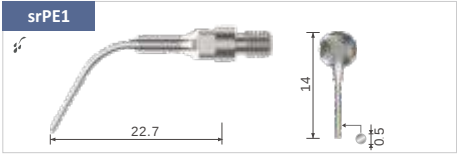

Perio. For rinsing and root cleaning in deep gum pockets.

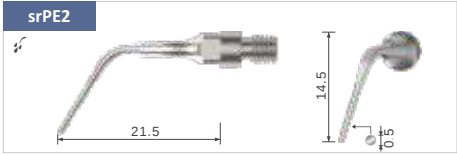

Perio. Curved to the right, for root cleaning in deep gum pockets in the posterior region.

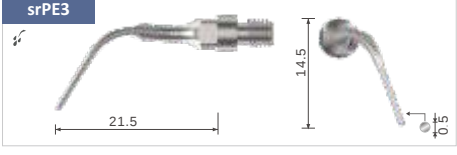

Perio. Curved to the left, for root cleaning in deep gum pockets in the posterior region.

IMPLANT CLEANING



Implant Cleaning. The ICP tip is used for sub- and supragingival cleaning of implants and restorations. The replaceable peek makes the action safe and gentle. The set contains 1 ICP, 1 ICP wrench and 2 spare peeks.



Replacement. ICPpeeks are autoclavable and reusable. Replacing the ICPpeek is simple using the wrench included in ICP SET.






Xpedent tips (suitable for use with KAVO® scaler machines).
The Xpedent brand belongs to a third party which has no connection with the original equipment manufacturer or its representatives.

Xpedent tips suitable for use with Sonicflex®


SCALING

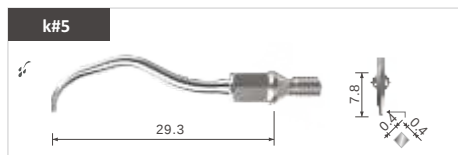

General. For calculus removal from all tooth surfaces.



Fine. Fine instrument for smaller areas and those which are difficult to access.

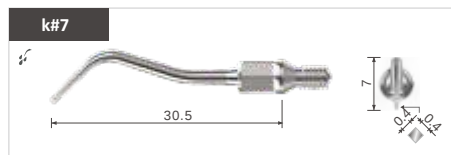

Subgingival. Instrument used for subgingival scaling.




Subgingival. For the removal of subgingival calculus.


PERIO



Perio. A long straight tip with rounded surfaces, the k#60 is used for gentle pocket cleaning during periodontal treatment. Used in the anterior region, it is effective on smooth surfaces and at cleaning pockets up to 9 mm deep.

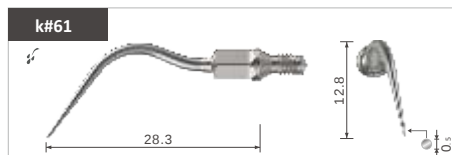


Perio Left. Curved to the left, for buccal access in the first and third quadrants to central incisor. Gentle but effective cleansing in pockets up to 9 mm deep, then cleansing of front teeth. Lingual access to premolars and molars in the second and fourth quadrants.

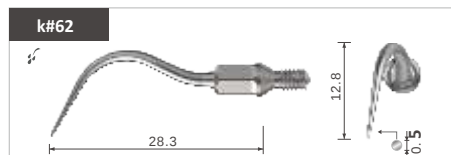


Perio Right. Curved to the right, for buccal access in the first and third quadrants to central incisor. Gentle but effective cleansing in pockets up to 9 mm deep, then cleansing of front teeth. Lingual access to premolars and molars in the second and fourth quadrants.



PerioSlim. The fast and minimally-invasive kPS instrument is ideal for calculus removal in deep periodontal pockets in all quadrants.

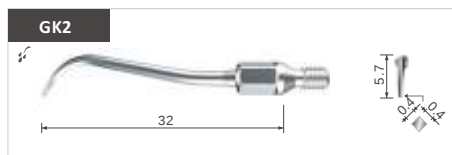


Used for supragingival, subgingival and interdental calculus removal.

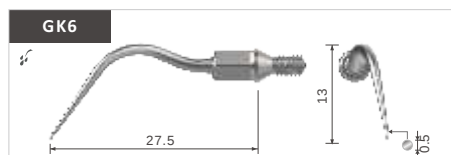

Curved to the left, used for supragingival, subgingival and interdental calculus removal.

GK7

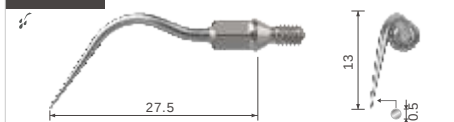

Curved to the right, used for supragingival, subgingival and interdental calculus removal.


EK8



FG-bur Holder. FG-bur holder used to expand the root canal and grind teeth.

ENDO


k#66

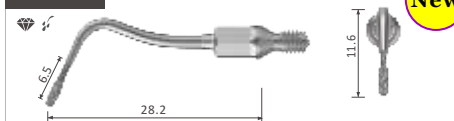

Access Shaping. Big bud shaped diamond coated tip for shaping the access cavity.

New


k#67

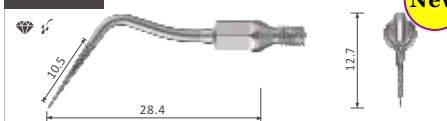

Canal Expansion. Diamond coated 125° conical tip. Used for widening canals.

New


k#68

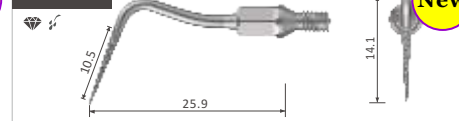

Canal Expansion. Diamond coated 112° conical tip. Used for widening canals.

New


k#69

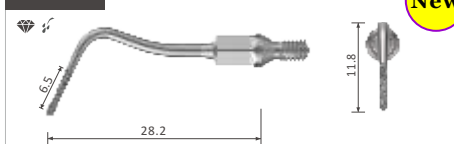

Canal Preparation. Small bud shaped diamond coated tip for removal of unwanted materials.

New

k#70

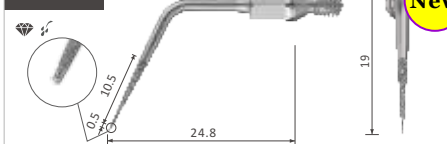


Long Canal Preparation. Long diamond coated 117° conical tip. Tip end not diamond coated for long and wide canals.

New

SURGERY


k100

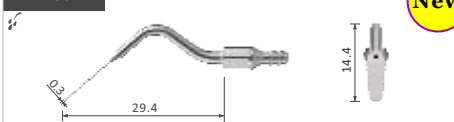

Minimally Invasive Surgery - Sagittal. Used for tooth detachment and extraction, fine bone cutting and splitting the alveolar crest. Easy to handle and gentle on soft tissues. Cutting depth of 10.7mm

New


k101

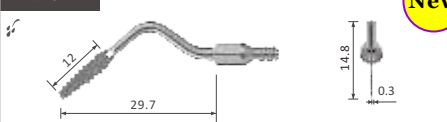

Minimally Invasive Surgery - Axial. Used for tooth detachment and extraction, fine bone cutting and splitting the alveolar crest. Easy to handle and gentle on soft tissues. Cutting depth of 10.7mm

New

k102

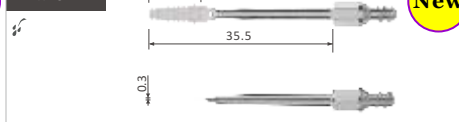

Minimally Invasive Surgery - Straight. Used for tooth detachment and extraction, fine bone cutting and splitting the alveolar crest. Easy to handle and gentle on soft tissues. Cutting depth of 10.7mm

New

IMPLANT CLEANING



Implant Cleaning. The ICP tip is used for sub- and supragingival cleaning of implants and restorations. The replaceable peek makes the action safe and gentle. The set contains 1 ICP, 1 ICP wrench and 2 spare peeks.




Replacement. ICPpeeks are autoclavable and reusable. Replacing the ICPpeek is simple using the wrench included in ICP SET.




Xpedent tips suitable for use with Sonicflex quick®


SCALING

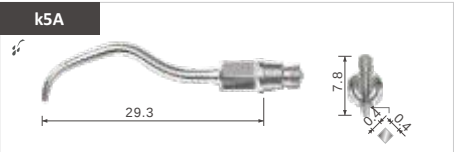


General. For calculus removal from all tooth surfaces.

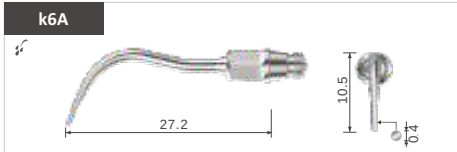

Fine. Fine instrument for smaller areas and those which are difficult to access.

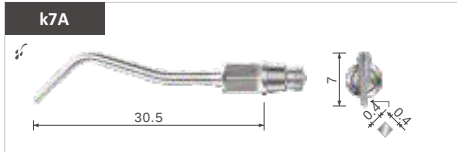


Subgingival. Instrument used for subgingival scaling.

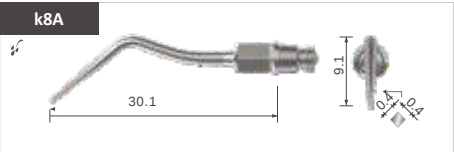

Subgingival. For the removal of subgingival calculus.

PERIO

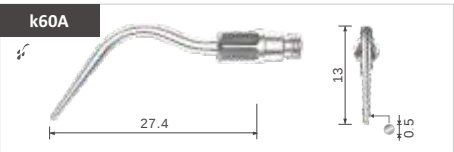

Perio. A long straight tip with rounded surfaces, the k60A is used for gentle pocket cleaning during periodontal treatment. Used in the anterior region, it is effective on smooth surfaces and at cleaning pockets up to 9 mm deep.



Perio LeftP. Curved to the left, for buccal access in the first and third quadrants to central incisor. Gentle but effective cleansing in pockets up to 9 mm deep, then cleansing of front teeth. Lingual access to premolars and molars in the second and fourth quadrants.




Perio Right. Curved to the right, for buccal access in the first and third quadrants to central incisor. Gentle but effective cleansing in pockets up to 9 mm deep, then cleansing of front teeth. Lingual access to premolars and molars in the second and fourth quadrants.




PerioSlim. The fast and minimally-invasive kPSA instrument is ideal for calculus removal in deep periodontal pockets in all quadrants.


ENDO

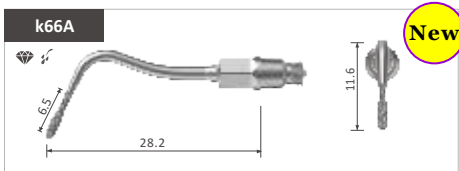


Access Shaping. Big bud shaped diamond coated tip for shaping the access cavity.

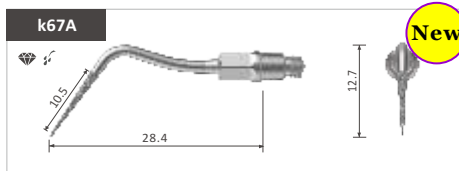


Canal Expansion. Diamond coated 125° conical tip. Used for widening canals.

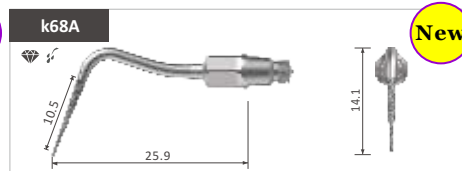

Canal Expansion. Diamond coated 112° conical tip. Used for widening canals.

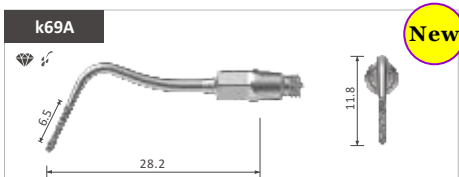

Canal Preparation. Small bud shaped diamond coated tip for removal of unwanted materials.

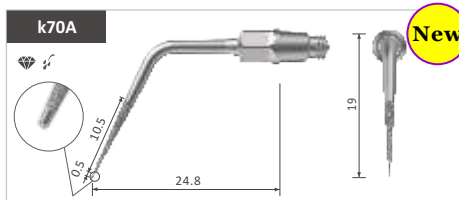

Long Canal Preparation. Long diamond coated 117° conical tip. Tip end not diamond coated for long and wide canals.

IMPLANT CLEANING



Implant Cleaning. The ICP tip is used for sub- and supragingival cleaning of implants and restorations. The replaceable peek makes the action safe and gentle. The set contains 1 ICP, 1 ICP wrench and 2 spare peeks.



Replacement. ICPpeeks are autoclavable and reusable. Replacing the ICPpeek is simple using the wrench included in ICP SET.


TW KAV Q



Torque wrench suitable for use with
Sonicflex®, Piezolutux®
& Sonosoft®, Sonicflex quick®

Xpedent tips suitable for use with Piezolux® & Sonosoft®


SCALING & PERIO

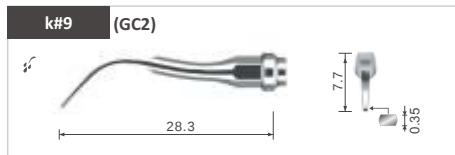

Sickle. Used to remove calculus and bacterial plaque from supragingival interdental calculus.



Perio. Used to remove calculus and bacterial plaque from supragingival interdental calculus.

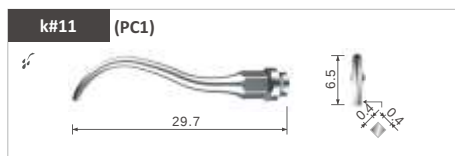


Universal. Used to remove subgingival calculus.

TW KAV Q




Torque wrench suitable for use with Sonicflex®, Piezolux® & Sonosoft®, Sonicflex quick®


Xpedent tips suitable for use with Piezoled® & Piezosoft®

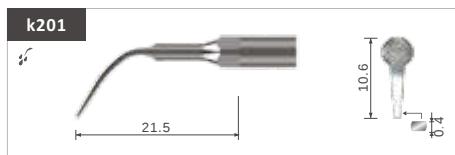

Universal. Accurate removal of solid supragingival calculus in all quadrants.

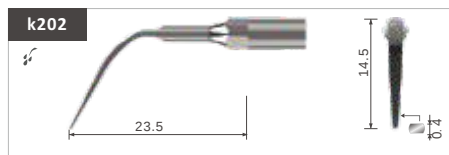

Perio. Well suited for the precise removal of stubborn subgingival and supragingival calculus and concretions in all quadrants.

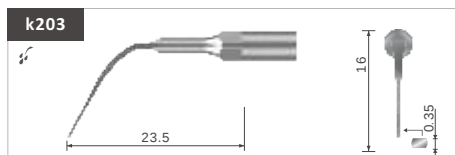

PerioSlim. Fast and minimally-invasive, ideal for calculus in deep periodontal pockets in all quadrants.



File Holder 120°. Use with Xpedent U files with a shank diameter of 0.57mm, this tip is usually used for anterior root canal cleaning.

TW KAV M



Torque wrench suitable for use with KAVO Piezoled® & Piezosofts®



Xpedent tips (suitable for use with AMDENT® scaler machines).
The Xpedent brand belongs to a third party which has no connection with the original equipment manufacturer or its representatives.

SCALING & PERIO




Universal. Tip for subgingival scaling, furcations, supragingival fine scaling and spot removal.



Perio. Primarily for lingual and buccal subgingival scaling and furcations. Also used for supragingival fine scaling and spot removal.



Heavy Calculus Removal. Can be used for lingual and buccal supragingival scaling.

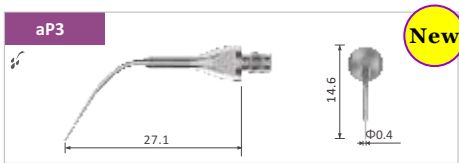

Perio. Used to remove subgingival calculus.

TW AMD M



















Suitable for use with
Amdent® tips








Xpedent tips (suitable for use with MECTRON® surgery machines).
The Xpedent brand belongs to a third party which has no connection with the original equipment manufacturer or its representatives.

BONE CUTTING

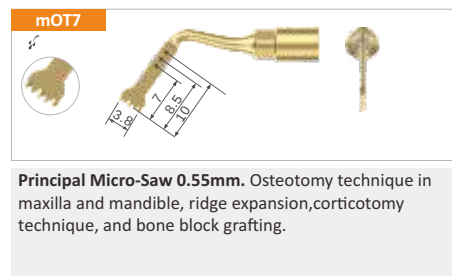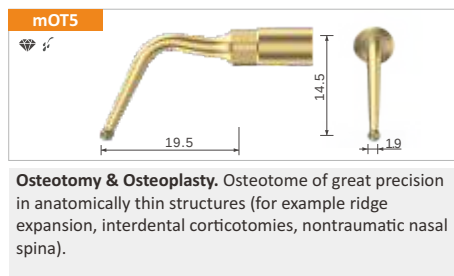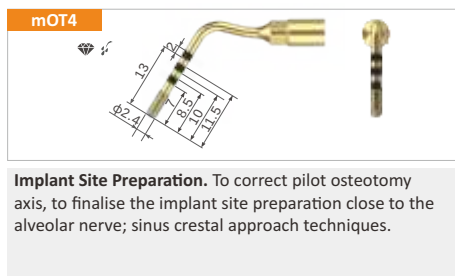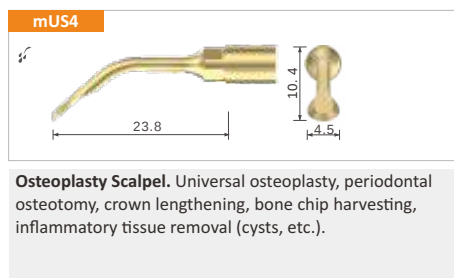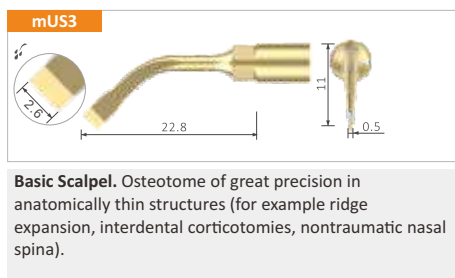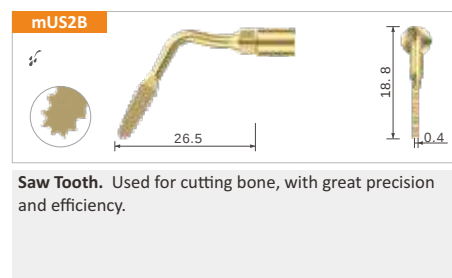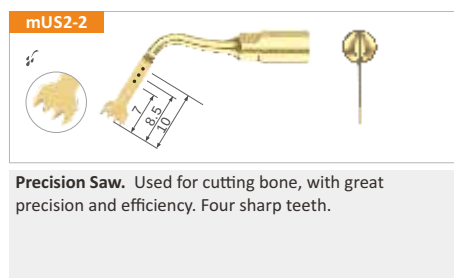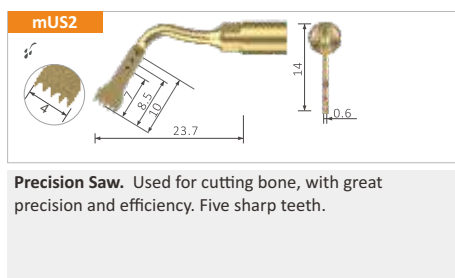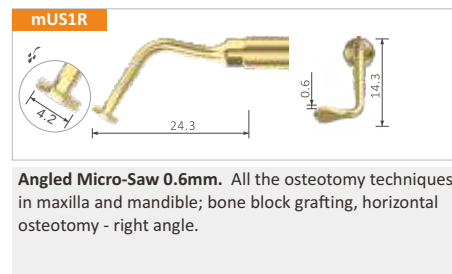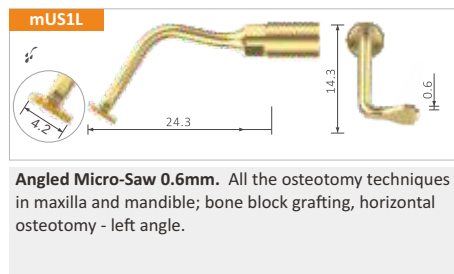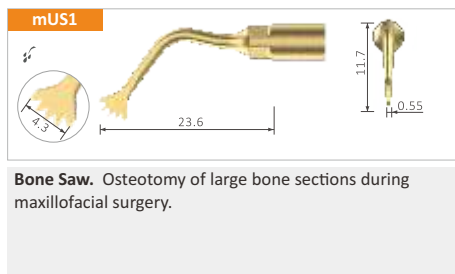

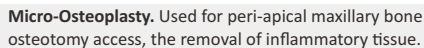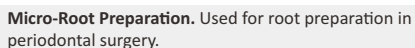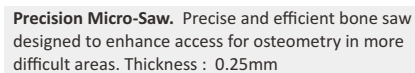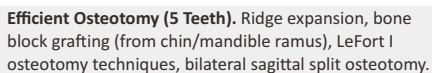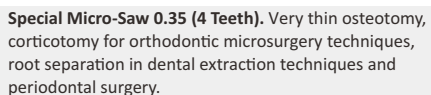


mUL1

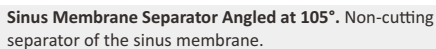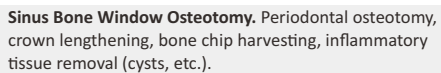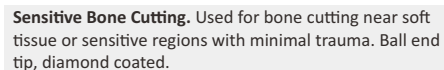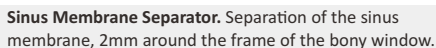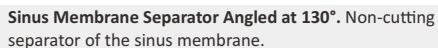


IMPLANTATION

mUI1

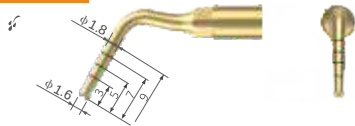


1.6 mm Enlargement. Used to enlarge implantation area avoiding soft tissue trauma.

mUI2

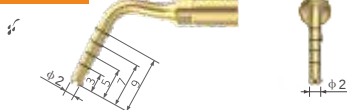


2.0 mm Enlargement. Used to enlarge implantation area avoiding soft tissue trauma.

mUI7

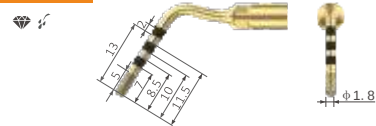


Fine Preparation. Used for final preparation of bone implant sites.

mUI8

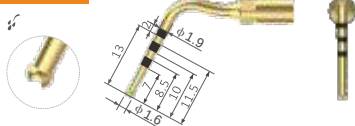

Implant Preparation. Diameter 1.6mm x 15mm. Used to enlarge implant site.

mUI9

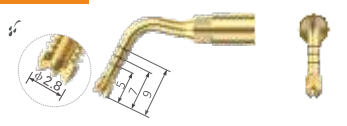


Implant Preparation. 2.8 mm Cutting head.

mUI9E




Deep Preparation. 2.8 mm Cutting head.

mIM1

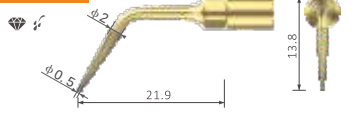


Diamond Bone Perforation. Conical surgery tip with diamond layer for bone perforation. Recommended for pilot osteotomy in extraction sites and to check the alignment of the axis of PIN IM1 preparations.

mIM2A

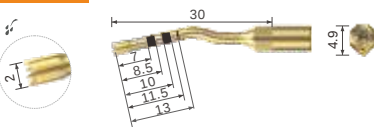


Bone Perforation. Preparation of the implant bed, especially bone perforation. Particularly effective for pilot osteotomy in the anterior and posterior (maxillary) region or to check the alignment of the axis of the preparations 2 PIN 2-2.4.

mIM2P

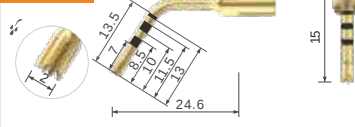


Bone Perforation. Pilot osteotomy in posterior region (in the maxilla).

mIM3A

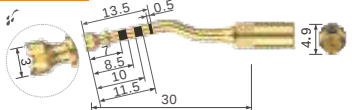


Bone Perforation. To enlarge or to finalise the implant site preparation; insert with double irrigation to avoid overheating (in the maxilla).

mIM3P

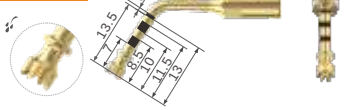


Bone Perforation. To enlarge or to finalise the implant site preparation; insert with double irrigation to avoid overheating (in the maxilla).

mIM4A

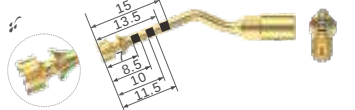


Bone Perforation. To finalise the implant site preparation; insert with double irrigation to avoid overheating (in the maxilla).

mIM4P

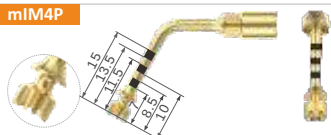


Bone Perforation. To finalise the implant site preparation; insert with double irrigation to avoid overheating (in the maxilla).

mIP2-3

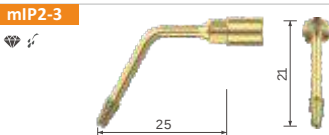

Pilot Implant Site Preparation. To optimise concentricity of implant site preparation between $\varnothing 2$ and $\varnothing 3$ mm.


mIP3-4

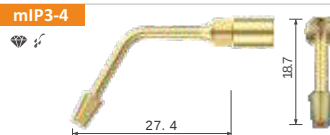

Micrometric Osteotomy. To optimise concentricity of implant site preparation between $\varnothing 3$ and $\varnothing 4$ mm.


RETRO

mUE1

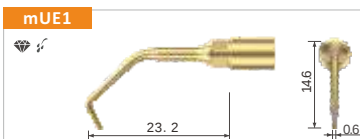


Efficient Canal Cleaning. Apical root debridement.

mUE2

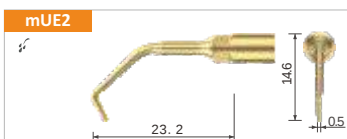


Gentle Canal Cleaning. Gentle apical root debridement.

mUE3

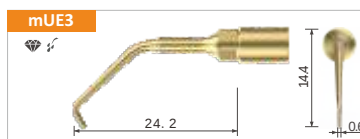

Efficient Canal Cleaning. Apical root debridement.


mUE4

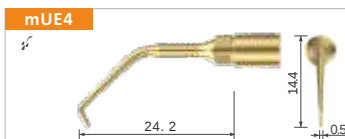

Gentle Canal Cleaning. Gentle apical root debridement.

EXTRACTION

mEX1

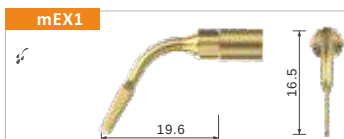

Extraction Scalpel. Root osteoplasty to cut off the ankylosis, root fraction techniques.

mEX2




Angled Extraction Scalpel. Root osteoplasty in the posterior regions.

mEX3



Angled Micro-Saw 0.6mm. Root osteoplasty in the posterior regions.

mXM-NINJA

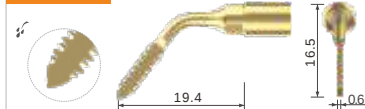

Tooth Sectioning. Saw-tooth tip particularly efficient for hemisections and root amputations.

mUC1



Extraction Scalpel. Root osteoplasty to cut off the ankylosis, root fraction techniques.


SCALING

mUP1



Gentle Scaling. Root scaling.

mUP2

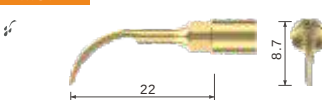


Periodontal Scaler. Scaling and inflammatory tissue removal, fractured root apex extraction.

mUP3




Angled Curette. Angled curette root scaling.

mUP4

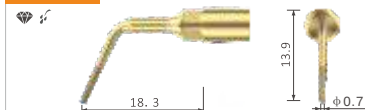

Scaling and Microdebridement. Root debridement and root planing during resective and regenerative periodontal surgery.

mUP5

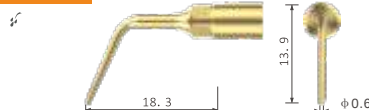


Root Surface Micro-Smoothing. Root planing.

mUP6



Micro-Root Preparation. Root preparation in periodontal surgery.

mUP7

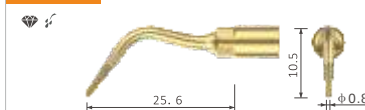

Micro-Osteoplasty. Interproximal osteoplasty and root planing.

IMPLANT CLEANING

mICP SET



Implant Cleaning. The ICP tip is used for sub- and supragingival cleaning of implants and restorations. The replaceable peek makes the action safe and gentle. The set contains 1 ICP, 1 ICP wrench and 2 spare peeks.

ICPpeek



Replacement. ICPpeeks are autoclavable and reusable. Replacing the ICPpeek is simple using the wrench included in ICP SET.

TW MEC M

New



Suitable for use
with Mectron® tips



Xpedent tips (suitable for use with EMS[®] surgery machines).
The Xpedent brand belongs to a third party which has no connection with the original equipment manufacturer or its representatives.

Xpedent tips (suitable for use with SATELEC[®] surgery machines).
The Xpedent brand belongs to a third party which has no connection with the original equipment manufacturer or its representatives.

BONE SURGERY



Cortical Bone Saw. Saw intended for cutting of cortical bone in surgical orthodontics, pre-implant surgery, and bone distraction.

E-BS1

BS1

V-BS1



Left Cortical Bone Saw. Left-oriented bone saw used for cutting the ramus cortical bone.

E-BS2L

BS2L

V-BS2L



Right Cortical Bone Saw. Right-oriented bone saw used for cutting the ramus cortical bone.

E-BS2R

BS2R

V-BS2R



Circular Scalpel. Circular scalpel at 130° for osteoplasty and harvesting of bone particles or chips.

E-BS4

BS4

V-BS4



Flat Scalpel. For fine osteotomies. All high-precision osteotomy techniques such as crest expansion, osteotomy line, distraction, approach of the vestibular bone flap in the presence of a thick cortex in sinus surgery.

E-BS5

BS5

V-BS5



Curved Scalpel. Effective for substantial osseous reshaping. It is used to perform osteoplasties and reshaping of the bone not assuring the tooth support. It can also be used to mark a prosthetic reference point on enamel, where osteoplasty will begin.

E-BS6

BS6

V-BS6



Special Micro-Saw. Osteotomy technique in maxilla and mandible, ridge expansion, corticotomy technique and bone block grafting.

E-BS6F

IMPLANTATION



Diamond Bone Perforation. Conical surgery tip with diamond layer for bone perforation. Recommended for pilot osteotomy in extraction sites and to check the alignment of the axis of PIN IM1 preparations.

XS-IM1

V-IM1



Bone Perforation. Preparation of the implant bed, especially bone perforation. Particularly effective for pilot osteotomy in the anterior and posterior (maxillary) region or to check the alignment of the axis of the preparations 2 PIN 2-2.4.

XS-IM2A

V-IM2A



Bone Perforation. To enlarge or to finalise the implant site preparation; insert with double irrigation to avoid overheating (in the maxilla).

XS-IM3A

V-IM3A

SINUS LIFTING



Sinus Bone Window Osteotomy. Diamond-coated tip for vestibular bone window cut and attenuation of sharp angles.

E-SL1

SL1

V-SL1



Micrometric Osteoplasty. Diamond-coated ball tip for smoothing the vestibular bone window and precise osteoplasty.

E-SL2

SL2

V-SL2



Sinus Membrane Separator. Non-cutting tip served for schneiderian membrane elevation on the window's edges.

E-SL3

SL3

V-SL3



Sinus Membrane Elevator. Non-cutting 90° angled spatula allows the sinus membrane to be lifted deep inside the sinus. Diameter of the spatula: 4mm

E-SL4

SL4

V-SL4



Non-Cutting Spatula. Non-cutting spatula used for schneiderian membrane elevation inside the sinus and for disengagement of anatomical structures.

E-SL5

SL5

V-SL5

EXTRACTION



Root Osteoplasty. Used for syndesmotomy and periradicular osteotomy.

E-LC1

LC1

V-LC1



Extraction Scalpel. Used for cutting ligaments with minimal trauma when pulling out firm and solid teeth.

E-LC1/90

LC1/90

V-LC1/90



Slim Scalpel. Slim tip, particularly slim, allows for access to the tight spaces between the root and the alveolar bone.

E-LC2

LC2

V-LC2



Angled Scalpel. Left-oriented at 45°, tip facilitates access to the posterior areas.

E-LC2L

LC2L

V-LC2L



Angled Scalpel. Right-oriented at 45°, tip facilitates access to the posterior areas.

E-LC2R

LC2R

V-LC2R



Root Amputator. Saw-tooth tip particularly efficient for hemisections and root amputations.

E-Pirahna

Pirahna

V-Pirahna

XPEDENT TORQUE WRENCHES & WRENCHES

TW UNI P



Universal
Plastic Torque Wrench for use
with EMS®, SATELEC®
and NSK® scaler tips

TW SIR P



Plastic Torque Wrench
suitable for use
with SIRONA® tips

TW KAV Q



Metal torque wrench suitable
for use with Kavo
Sonicflex®,
Piezolux® & Sonosoft®
and Sonicflex quick®

TW KAV M



Metal torque wrench suitable
for use with Kavo
Piezoled® & Piezosoftware®

TW MEC M



Metal torque wrench suitable
for use with
all surgical Mectron® tips

TW AMD M



Metal torque wrench suitable
for use with
all Amdent® tips

New

W ENDO



Endo Wrench

ICP WREN



ICP Wrench

KICP WREN

New



ICP Wrench

SIRONA® is a registered trade mark of Sirona Dental Systems GmbH
EMS® is a registered trade mark of Electro Medical Systems SA
SATELEC® is a registered trade mark of ACTEON Group
Mectron® is a registered trade mark of Mectron s.p.a.
Amdent® is a registered trade mark of LM Instruments Oy
Kavo® is a registered trade mark of KAVO Dental Corporation
NSK® is a registered trade mark of Nakanishi INC. Corporation

Sonicflex®, Piezolux® & Sonosoft®, Piezoled® & Piezosoftware®, Sonicflex quick® are registered trade marks of the company KAVO DENTAL GMBH

U-FILES (32mm)

UF-SET/32	UF-15/32	UF-20/32	UF-25/32	UF-30/32
				
U-Files (32mm) This set contains 1 piece each of UF-15/32, UF-20/32, UF-25/32 and UF-30/32.	ø0.15mm 10pieces/box	ø0.2mm 10pieces/box	ø0.25mm 10pieces/box	ø0.3mm 10pieces/box

Attention: Our U-Files work with the following file holders:
E1, ES1, EN1, ED1, E2, ES2, EN2, ED2 and k222.

THREAD IDENTIFICATION

TC-1 	EMS® - Surgery	TC-4 	MECTRON® - Surgery KAVO PiezoLED® / PIEZOsoft®
TC-2 	EMS® MECTRON® KAVO®-PIEZOlux® / SONOsoft® SATELEC® NSK®	TC-5 	SATELEC® - Piezo1 SATELEC® - Piezo2
TC-3 	KAVO® - SONICflex® NSK® - Air	TC-6 	SIRONA® AMDENT®

SIRONA® is a registered trade mark of Sirona Dental Systems GmbH
EMS® is a registered trade mark of Electro Medical Systems SA
SATELEC® is a registered trade mark of ACTEON Group
Mectron® is a registered trade mark of Mectron s.p.a.
Amdent® is a registered trade mark of LM Instruments Oy
Kavo® is a registered trade mark of KAVO Dental Corporation
NSK® is a registered trade mark of Nakanishi INC. Corporation
Sonicflex®, Piezolux® & Sonosoft®, Piezoled® & Piezosoft®, Sonicflex quick® are registered trade marks of the company KAVO DENTAL GMBH

XPEDENT SCALER HANDPIECES

HP YE



NO-LED Handpiece suitable
for use with EMS® scaler machines

HP YS



NO-LED Handpiece suitable
for use with SATELEC® scaler machines

HP YE7 LED



LED Handpiece suitable
for use with EMS® scaler machines

HP YS8 LED



LED Handpiece suitable
for use with SATELEC® scaler machines

Xpedent handpieces are manufactured from aluminium, they are therefore able to withstand autoclaving without the exterior cracking problems often associated with plastic handpieces. They are suitable for use with EMS® and SATELEC® machines. LED handpieces have a different connection from non-LED.

AUTOCLAVABLE TIP HOLDERS

TH WC



Tip holder with cover

TH



Tip holder

ABOX purple

New




Tip holder with cover

XPEDENT BUILT-IN SCALER MACHINES



HP YE7 LED

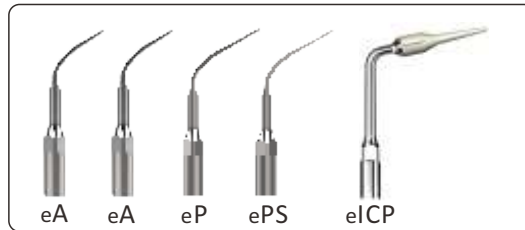

5 Tips



TW UNI P



SHEL-BE YE

SHEL-BE YS

Built-in Scaler: For installation into dental chair.

TIP WALLETS

DEMO1



DEMO1S



DEMO3



DEMO4



The tip wallets are made from black anodised aluminium. They are designed to show our tip range to best advantage. The standard wallets DEMO1 and DEMO1S contains our recommended range. Wallets can be customised for ODM or OEM use. They are recommended for both use at dental shows and by travelling sales staff.



depuis 2011

L'EXCELLENCE AU SERVICE DES CABINETS DENTAIRES

CONTACTEZ-NOUS OU COMMANDEZ SUR
NOTRE SITE INTERNET POUR BÉNÉFICIER
DE CES OFFRES EXCLUSIVES.

PRODENTALPES

 06 01 21 16 35

 contact@prodentalpes.fr

 ZI DES DRAGIEZ, 30 RUE DE GRÉBELIN,
74800 LA ROCHE-SUR-FORON



10 000 PRODUITS TOUTES MARQUES

Totalement indépendant des laboratoires, nous sommes à votre écoute
et nous vous conseillons de façon totalement objective.

Visitez le site prodentalpes.fr



Découvrez nos Promotions & Formations